

NEW

# Introducing our new SCFvh 4032 combination



No Frost

## NoFrost technology

With professional laboratory appliances from Liebherr, defrosting is no longer necessary. Any moisture that occurs is bound by the fan-assisted cooling, and needs-based defrosting phases drive even the last moisture particles out of the interior. For convenient, safe, and ice-free storage.

## Unique features



Ultra Efficient

Energy-efficient refrigerators and freezers help with ever increasing electricity prices. Liebherr fully utilizes the savings potential of every component – from the compressor to the seal, to the lighting and electronics. The ultra-efficient appliances in the Perfection series for example demonstrates how this has a positive effect on your overall operating costs.



Laboratory and pharmacy refrigerators from Liebherr meet all the necessary equipment and performance criteria to fully comply with DIN 13277. DIN 13277 defines the requirements and test specifications for refrigerators and freezers used in laboratories and for storing medicinal products.



**+2 °C SafetyDevice**  
If the temperature in the appliance drops below +2 °C due to a malfunction, the +2 °C SafetyDevice function takes control, which keeps the temperature stable and activates an alarm. Sensitive substances are protected from freezing and you are notified immediately.



**Mechanical Lock**  
Liebherr professional appliances are equipped with a mechanical lock and two keys. This ensures that unauthorised persons cannot access your valuable and sensitive stored goods.




**Environmentally friendly refrigerants**  
Liebherr appliances are optimized for maximum energy efficiency and low operating costs. We only use the natural and environmentally-friendly refrigerants R 600a and R 290. When used with highly efficient compressors, these are extraordinarily efficient and, with their low global warming potential, also future-proof.

Learn more at: [home.liebherr.com](http://home.liebherr.com)

**LIEBHERR**

# SCFvh 4032 Combination

13.7 cubic feet (3.9 Freezer/ 9.9 Refrigerator)

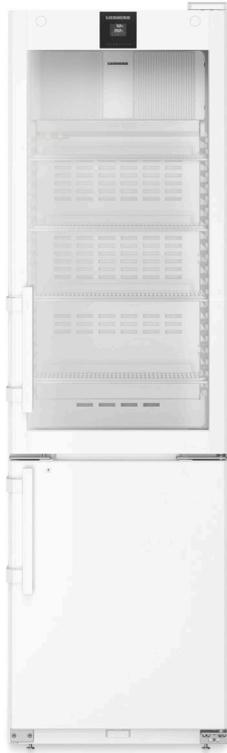

Total gross / net capacity l (cu. ft.) refrigerator	9.9 / 8.7
Total gross / net capacity l (cu. ft.) freezer	3.9 / 2.8
Exterior dimensions in mm (w / d / h) in inches (w / d / h)	597 / 654 / 2044 23 9/16   25 3/4   80 1/2
Ambient temperature range °F	*+50°F to +95°F
Refrigerant	R 600a
Sound rating**	48 dB
Connection rating	2,0 A
Frecuency / Voltage	60 Hz / 115-127V~
Cooling System / Defrost	Dynamic / Automatic
Temperature range in °C (°F) cooling	*+3°C to +16°C (+38°F to +60°F)
Temperature range in °C (°F) freezing	*+9°C to -30°C (+15°F to -22°F)
Door / lid material	Insulated glass door, steel
Interface	Wifi / LAN
Temperature display	External digital
Interior Light	LED light column, left
Adjustable shelves refrigerator / freezer	4/5
Shelf material	Grid shelves plastics coated
Shelf load weight in kg (lbs) refrigerator / freezer	45 / 40
Handle	Antimicrobial handle with opening mechanism
Lock	Mechanical
Self-closing door	Yes
Door hinges	Right reversible
Insulation in mm (in)	45 / 26
Gross / net weight in kg (lbs)	106 (233) / 100 (220)

